



Kit RAL Quick TB

Kinyoun method

Staining with Kit RAL Quick TB allows to detect mycobacteria, also called A.F.B (Acid-Fast Bacilli) in a fast and safe way.

Kit Quick TB is a 2-step staining technique:
RAL Carbolic Fuchsin stains mycobacteria pink.
Armand Solution stains the background of the preparation blue.
The use of Carbolic Fuchsin avoids to heat the slide and breathe in hazardous fumes.

Easy to use, Kit RAL Quick TB contains 2 pre-calibrated solutions to insure standardized results.
The 2 bottles are equipped with a dropper pipette to dispense the stain directly on the slide.


KEY FEATURES

Streamline your workflow


Work safely

Save time


Give reliable results

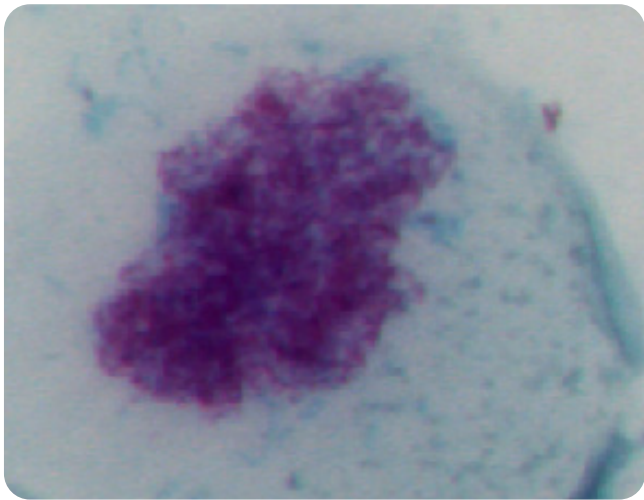


Objective x 40

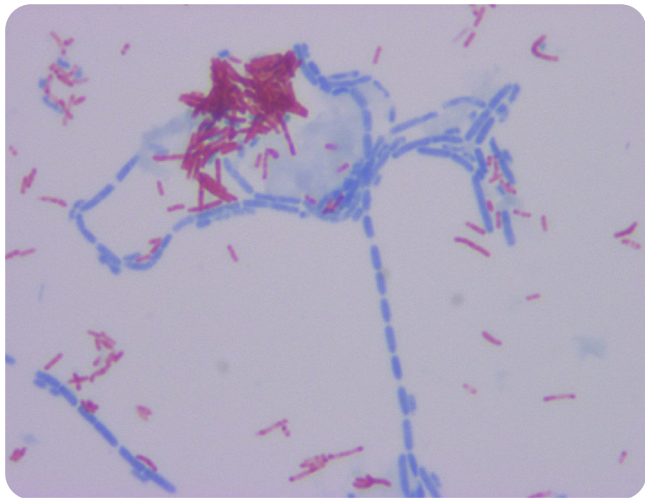

XXXX

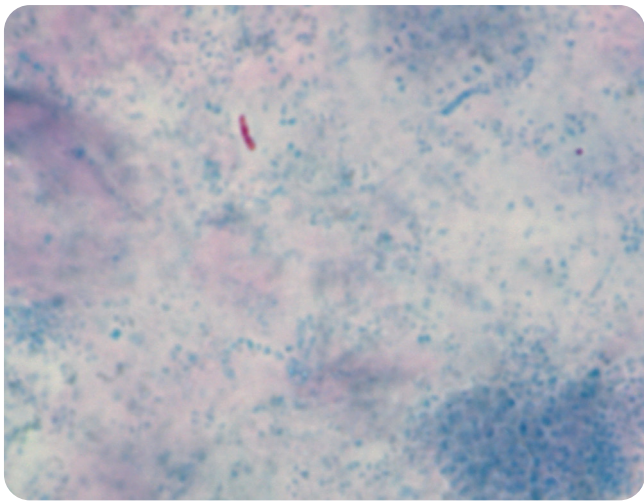

XXXXXX

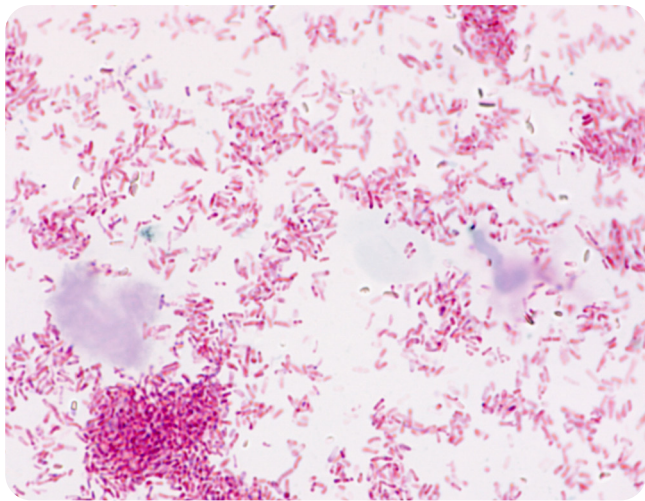

XX

Kit RAL Quick TB

AVAILABLE KIT



Kit RAL Quick TB REF.: 361560-0000

TECHNICAL SPECIFICATIONS

APPLICATION AREAS

- Bacteriology

NUMBER OF TESTS

- 200 slides

PROTOCOL DURATION

- 6 minutes

NUMBER OF BOTTLES

- 2 bottles

PRODUCT DESCRIPTION

- Carbolic fuchsin RAL: Clear dark red solution
REF.: 365240-0125
- Armand solution: Clear dark blue solution
REF.: 3601002A125

REFILLS

- Carbolic fuchsin RAL:
REF.: 365240-1000- 1x1L
- Armand solution:
REF.: 360100-1000- 1x1L